

# SYNCHRONIZED RF & HIFEM: FAT HISTOLOGY & SCANNING ELECTRON MICROSCOPY STUDY

SIMULTANEOUS APPLICATION OF HIGH-INTENSITY FOCUSED ELECTROMAGNETIC FIELD AND SYNCHRONIZED RADIOFREQUENCY FOR FAT DISRUPTION: HISTOLOGICAL AND ELECTRON MICROSCOPY PORCINE MODEL STUDY

Robert A. Weiss MD, FAAD<sup>1</sup>, MVDr. Jan Bernardy Ph.D.<sup>2</sup>,

Prof. MVDr. Frantisek Tichy, CSc<sup>3</sup>




1. Maryland Laser Skin & Vein Institute, Hunt Valley, MD, USA; 2. Veterinary Research Institute, Brno, CZ

3. Department of Anatomy and Histology, University of Veterinary and Pharmaceutical Sciences Brno, CZ

Published in Dermatologic Surgery journal, June 2021, DOI: 10.1097/DSS.0000000000003091

## HIGHLIGHTS

- Both histology and scanning electron microscopy showed damaged adipocytes post-treatment due to apoptosis and lipolysis.
- Adipocyte size was decreased by 31.1% at 2 weeks post-treatment.
- The temperature in fat tissue was maintained just below 45°C for the entire treatment.
- No necrosis was seen in the tissue.

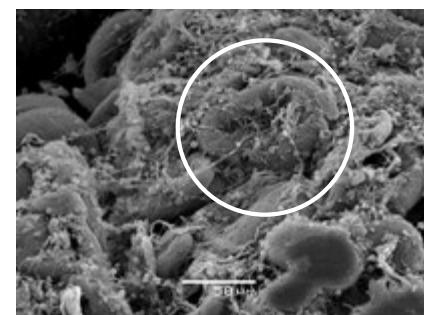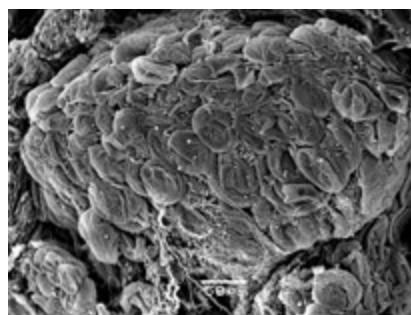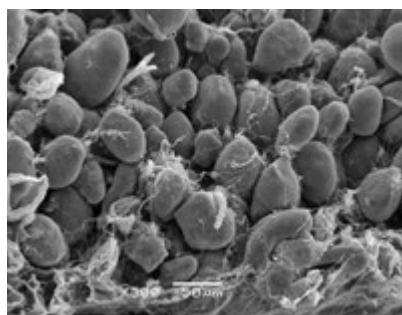



Healthy fat cells with well-defined shape at the baseline (left); shrunk adipocytes with noticeable membrane ruptures occurred at 4 days (center); disrupted adipocytes with extrusion of lipid droplets at two weeks (right)

## STUDY DESIGN

- 7 Large White pigs (approximately 6 months old).
- All animals received three 30-minute treatments applied to abdomen.
- Biopsies of fat tissue (including control specimens) were collected at baseline, 4 days, 2 weeks, 1 month and 2 months post treatment.
- Histologies were evaluated under the light and scanning electron microscope.

## CONCLUSION

- The procedure elevates the **temperature** in subcutaneous fat to levels **necessary for apoptosis induction**.
- **Efficacy of the procedure for adipocytes deletion** was documented in **252** analyzed tissue slices.
- Mild inflammatory response was present to promote the **apoptotic death cells removal**.
- **No burns, necrosis or other serious adverse events** were observed.

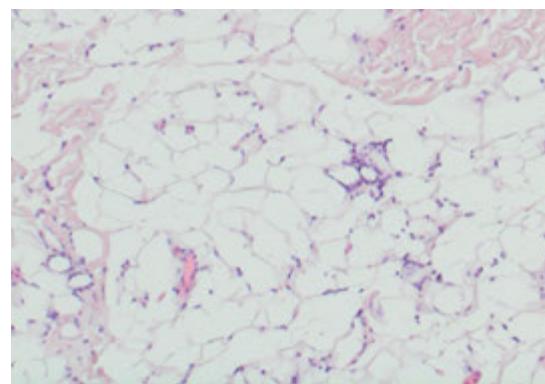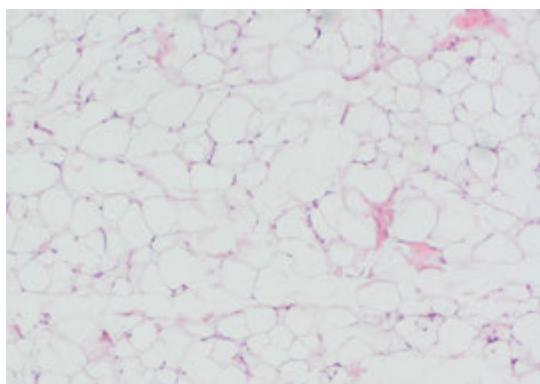

Baseline histology (left) showed adipocytes without any damage. At 2 weeks (right), flattened adipocytes with delaminated membranes are seen along with immune cells clearing the damaged tissue.